

VeraPlug™

Sterile Preloaded
Punctal Occlusion System

Description

The VeraPlug™ punctal occluder is designed to provide reduction or elimination of tear drainage through the inferior or superior puncta, thus maintaining lubricating tears on the surface of the eye. Each VeraPlug™ punctal occluder is molded from medical grade silicone. The VeraPlug™ is available in four sizes (small, medium, large and x-large) and is packaged sterile, two per box. Each occluder is sterile, preloaded on an inserter.

Indications for Use

The VeraPlug™ is for use in patients with dry eye syndromes.

Contraindications

Contraindications include, but are not limited to, eye infections, sensitivity or allergies to the occluder material and/or materials used in the manufacture of the device, blockage/infection of the lacrimal systems, inflammation of the eyelid, and epiphora.

Precautions

The VeraPlug™ may enhance the effect of some ocular medications in the eye. Depending on the type of medication being used the dose may need to be altered accordingly. If the patient experiences irritation, infection or epiphora after the insertion of the VeraPlug™, the occluder should be removed.

Potential Adverse Events

The following complications may occur:

- Epiphora
- Foreign body sensation
- Plug dislodgement or migration possibly requiring surgical intervention
- Pyogenic granuloma
- Infection of the lacrimal system
- Washout
- Punctal erosion

Product Features

Each box contains two individually packaged sterile VeraPlug™ punctal occluders preloaded on inserters for single use only. The VeraPlug™ is manufactured from implant grade silicone.

Proper Sizing

The VeraPlug™ punctal occluder is available in four sizes: small, medium, large and x-large. To determine the proper size, begin with the smallest VeraPlug™ gauge (0.5mm) and insert the head of the gauge into the punctal opening. The gauge should fit snugly with a small amount of resistance. If the gauge enters the punctal opening with no resistance, the next larger size gauge should be tried in the same manner. Repeat this process until the proper size VeraPlug™ is determined based on the sizing chart shown.



PUNCTAL OPENING	BEST FITTING GAUGE	PROPER VERAPLUG™ SIZE
0.4mm to 0.6mm	0.5mm	VSP-1001 Small
0.6mm to 0.7mm	0.6mm	VSP-1002 Medium
0.7mm to 0.8mm	0.7mm	VSP-1003 Large
0.8mm to 1.0mm	0.8mm	VSP-1004 X-Large

Prior to Insertion

Patients with intermittent tearing should receive probing and irrigation with a sterile irrigation solution to rule out pre-existing canalicular obstruction.

Insertion

- 1 Anesthetize the area of the punctum with a topical anesthetic placed in the conjunctival sac.
- 2 Apply a drop of saline solution or artificial tears onto the VeraPlug™ to help ease insertion.
- 3 Position the insertion instrument by placing the forefinger on the release button of the inserter and placing the occluder end of the insertion instrument over the patient's (superior or inferior) punctum.
- 4 Vertically insert the VeraPlug™ by positioning the occluder into the punctum until the cap is flush with the punctal opening. **FIGURE A**
- 5 When the occluder is properly seated, depress the release button and withdraw the insertion instrument. **FIGURE B**
- 6 Verify that the occluder is properly placed by confirming that the cap is flush with the punctal opening. **FIGURE C** After insertion, monitor the placement and integrity of the occluder to determine if/when the occluder may need to be replaced.

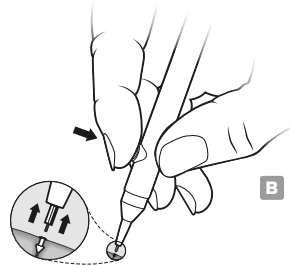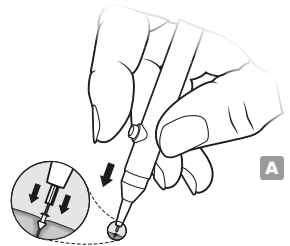

Removal

Should removal be indicated, grasp the vertical shaft of the occluder underneath the exposed cap with sterile forceps. Gently pull upward until the plug is removed. **FIGURE D**

Sterilization

VeraPlug™ sterile preloaded punctal occluders are offered in individual trays, two per box. The date of expiration should be confirmed prior to use. If the expiration date has lapsed the occluder should be discarded.

Storage

Store at room temperature.

Warnings

The VeraPlug™ punctal occluder is intended for single use. Do not reuse. If the sterile packaging is damaged or opened sterility is not guaranteed and the VeraPlug™ should be discarded. U.S. federal law restricts the sale of this device by or on the order of a physician.

Documentation

Each box contains Instructions for Use and two labels for your ease and convenience.

Authorized Representative

EC REP Emergo Europe
Prinsessegracht 20, 2514 AP
The Hague, The Netherlands

